


Microdialysis Fraction Collector, FC-90

The new FC-90 Fraction Collector saves both lab space and your time when collecting microdialysis samples from many animals, up to 4 at once. The FC-90's quiet operation and narrow (7" wide) footprint are ideal for small spaces.

The FC-90 is easy to program and controlled by Android™ devices via Bluetooth. The FC-90 has a robust design for microdialysis.

The newly designed Double Needles prevent wear and contamination. The outer stainless steel needle pierces seals and septa, protecting the inner PEEK needle from wear. The inert PEEK needle delivers the sample to the bottom of wells without exposing sample fluids to metal components. Refrigerated to 6°C without condensation issues on the samples. The eliminated cooling plate allows using various types of plates.

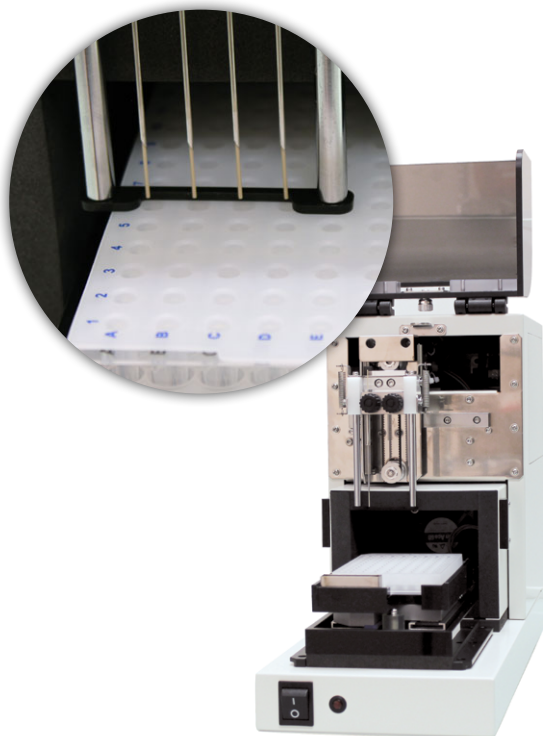

FC-90 Specifications	
Channels	1, 2, or 4 simultaneous samples
Sample Capacity	96 well plate
Cooling	6°C to room temperature, Peltier cooling
Dimensions	7" (W) x 18" (D) x 13" (H) Clearance of 5" behind the unit recommended..
Sample Delivery	Double-Needle: outer stainless steel needle pierces seals and septa, inner PEEK needle delivers the sample to the vial or well
Instrument Control	Programmable by Android devices
I/O	Bluetooth